

Bed Head Units

Description

The wall mounted bed head unit is primarily designed for one bed but can be modified to fit several beds in a row and is custom made.

It is made of extruded aluminium profiles that are epoxy polyester painted in various RAL colours. The design uses a minimal profile yet provides a large interior space, which gives the potential for varied configurations making the system functional, ergonomic and easy to maintain.

The ARDYS system provides for the installation of electric and gas supplies in isolated chambers that meet all of the relevant international standards.

The aluminium profiles are resistant to scratches and cleaning products. Clip in polycarbonate diffusers and clip in aluminium end caps give the system a clean well fitted appearance and avoids the need for visible screws.



Features

- High lighting efficiency giving long-life cycle and low power consumption.
- Able to accommodate various Medical gas outlets.
- Able to accommodate electrical sockets and earth connections.

Options Available




- Ambient Light
- Reading Light
- Night Light
- Gas piping following any Standard (DIN, BS, AFNOR...)
- High current chamber
- Low current chamber
- Access to Terminal bars

Installation

The bed head unit is wall mounted by means of included steel brackets.

Specifications:

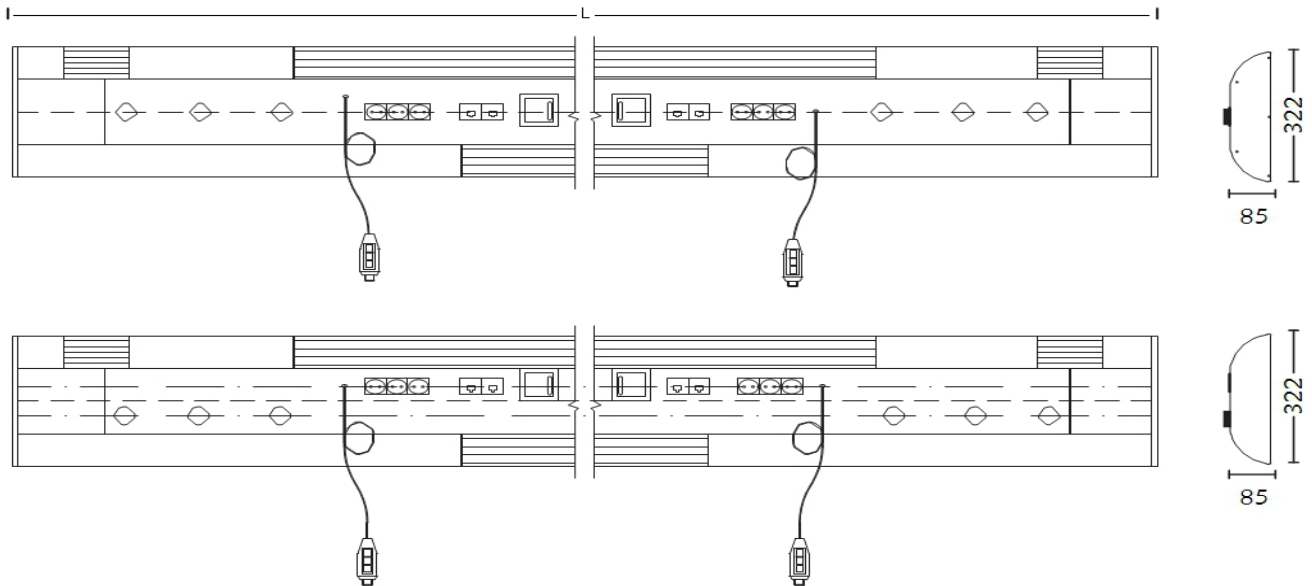

Electrical sockets: As required

Earth connections: As required

Lighting options:

- Ambient Light: One/two T5 bulbs located on the upper side of the BHU
- Reading Light: One/two T5 bulbs located on the lower side of the BHU
- Night Light: One T5 bulb located on the upper side of the BHU

Gas outlets: Located on the front side of the BHU. These cut-outs and holders can be modified during the installation on site.

Other options: Nurse call systems, phone jacks, microphone, speakers.



www.starkstrom.com

256 Field End Road, Eastcote, Ruislip, Middlesex, HA4 9UW, UK

Tel +44 (0) 20 8868 3732

Fax +44 (0) 20 8868 3736

