



STARKSTROM Ltd


wTCP Membrane Theatre Control Panel

The wTCP uses the very latest technology to provide hospitals with the flexibility and functionality not available from previous generations of Theatre Control Panels.

Advantages of the wTCP include:

- Modern contemporary styling.
- Wipe-clean – the panel facia is a continuous sheet of “Lexan” with no edges or grooves. Cleaning of the panel is therefore very easy.
- No more free-issue equipment on the panel; the wTCP easily interfaces with other systems within the theatre such as medical gas alarms, theatre lights, IPS and UPS.

The picture opposite shows a typical arrangement for alarms and controls on the wTCP. If, for example, there is a problem with the oxygen supply, an audible alarm will sound and the green “O₂ OK” will be replaced with a red “Fault” LED. The programming of when the system re-alarms is user-definable.

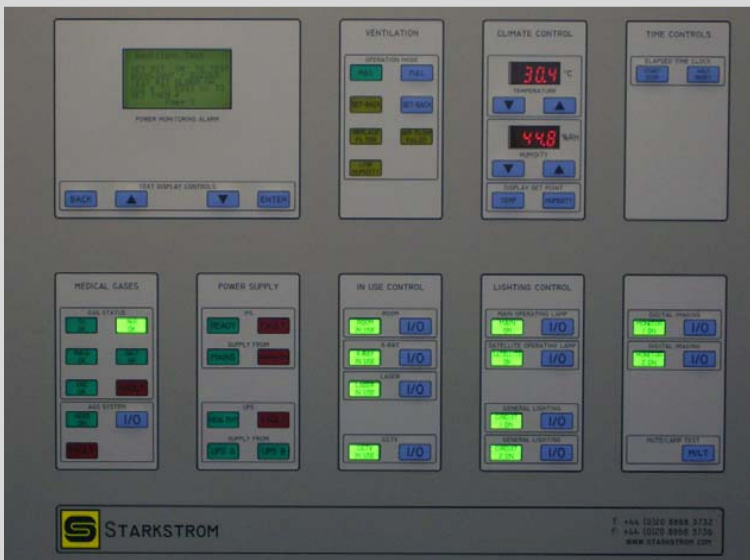

Typical dimensions of the wTCP are 897mmH x 597mmW x 120mmD. The panel is usually supplied with an analogue time-of-day and digital time elapsed clock. Other combinations of clocks can be offered together with matching X-Ray viewers and PACS monitors.



www.starkstrom.com

256 Field End Road, Eastcote, Ruislip, Middlesex, HA4 9UW, UK
Tel +44 (0) 20 8868 3732 Fax +44 (0) 20 8868 3736